




ARCIOLANE
2 viscosities



1300 cts
5500 cts

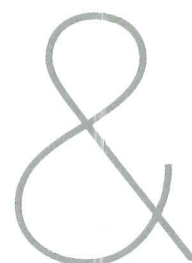

ARCIOLANE
2 presentations



Vial
Syringe



THE ASSURANCE
OF FRACTIONATED AND PURIFIED OIL



THE ASSURANCE
OF A REINFORCED SAFETY

The fractionation process removes the impurities and the low molecular weight polymers that clinical studies identified as particularly toxic.


Guarantees a quality conforming to the most demanding standards

	FDA requirements*	ARCIOLANE raw material
Low Molecular Weights below D9 or 666 g.mole ⁻¹	< 200 ppm	< 150 ppm
Resistivity	> 1. 10 ^{e14} ohm. cm	>1.10 ^{e14} ohm. cm
Volatility	<1 %	<0,1%

*based on FDA approved product characteristics

Guarantees a better tolerance

	Raw oil	ARCIOLANE
Low molecular weights content D4 to D9	148 ppm	21 ppm
Polydispersity index	2,04	1,86

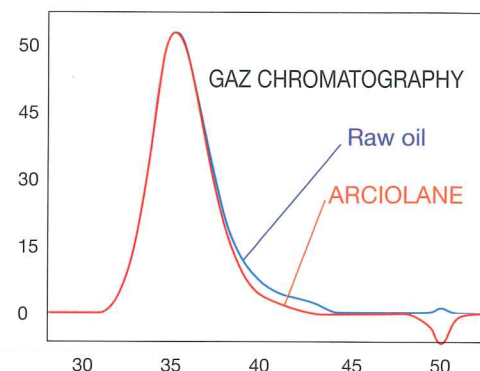

NO CHLOROBUTYL NO LATEX

PURE SILICONE STOPPER

NO EXTRACTABLE PHENYL COMPOUNDS

PHARMACOPOEIA COMPLIANT MOULDED GLASS

SHOCK-PROOF

STAINLESS STEEL LUER-LOCK

NO RISK OF DISCONNECTING



A UNIQUE CHOICE FOR ACCESSORIES





STERILE - SINGLE USE

ARC connection line for all vitrectomy equipment

Quick and easy :

- Provided with a 20 g blunt cannula : • 1 place the stopper into the syringe
• 2 insert the cap enclosing the glass flange • 3 connect the line to the equipment

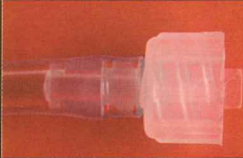

Reference	Type of Connector		Compatible with
ARC 510	Fixed Male / Female Luer-lock		Associate * - Harmony* *(Dorc) OpenPhaco (Corneal) Premiere (Storz) - Sto (Alcon) Penthaxis (Ioltech) - Antarès (Opticon)
ARC 520	rotative		Accurus (Alcon)
ARC 530	rotative		Millenium (B&L)
ARC 540	steady		Orbit (Oertli)

*Requires A 571 adaptor

APIS ARCIOLANE syringe
connection line for all
vitrectomy equipment



REUSABLE, autoclave
STERILIZABLE

Compatible with		
APIS alone	Premiere (Storz) - Sto (Alcon) - Penthaxis (Ioltech) - Antarès (Opticon)	
APIS + A 571	OpenPhaco (Corneal) - Associate (Dorc) Harmony (Dorc)	
APIS + A 572	Accurus (Alcon)	
APIS + A 573	Millennium (B&L)	
APIS + A 574	Orbit (Oertli)	

200A ANG-04 2011-04

mis manual ARCIOLANE syringe injection holder :
Easy, cost effective full control during the injection, can be re sterilized



ARCAD
O P H T A®

Find updated bibliographies on the website www.arcadophta.com

Arcadophta 11, rue Antoine Ricord
31100 TOULOUSE - FRANCE
Tél. : +33 5 61 40 52 35
Fax : +33 5 61 40 84 66
E-mail : info@arcadophta.com
Skype : arcad4