


## Specula


Adjustable or Non-adjustable Infant and Adult

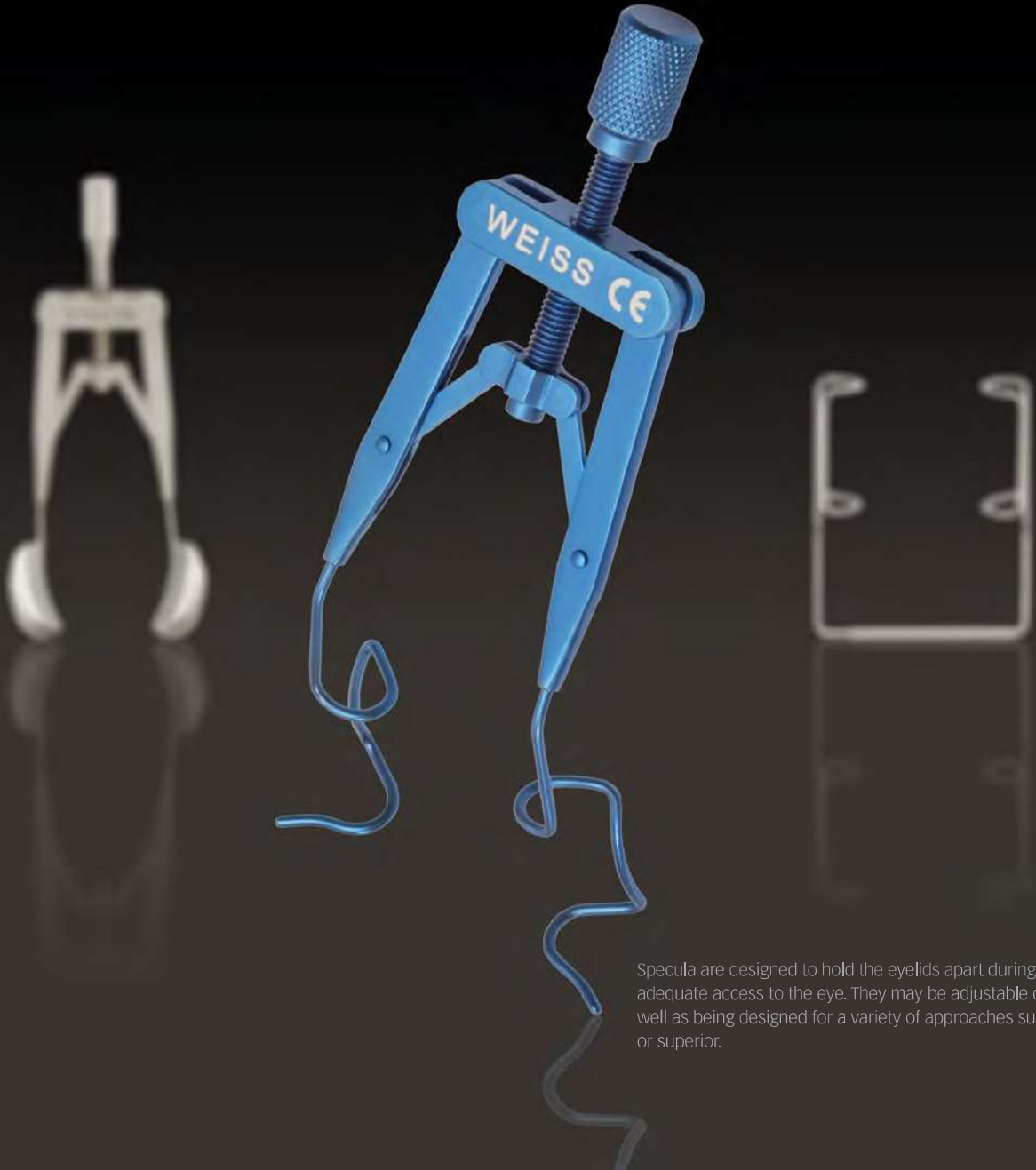


Specula are designed to hold the eyelids apart during surgery allowing adequate access to the eye. They may be adjustable or non-adjustable as well as being designed for a variety of approaches such as temporal or superior.

## Specula anatomy


### NON-ADJUSTABLE

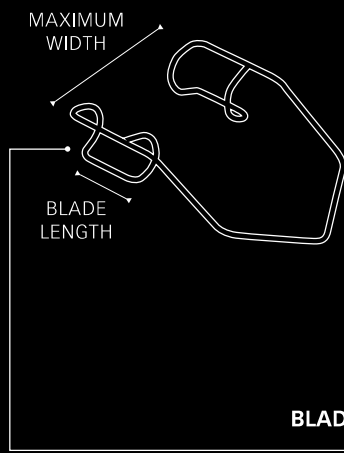

### ADJUSTABLE

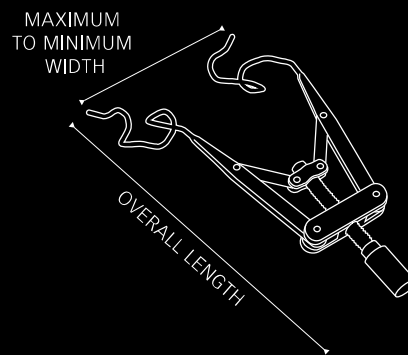

**Open:** allows good access for instruments  
e.g. Phaco handpiece that might otherwise rest on the top bar of a closed blade making surgery difficult

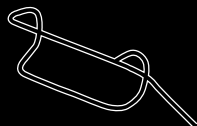

**Closed:** good hold on eyelids for procedures where the top bar does not cause a problem  
e.g. trabeculectomy or temporal incision phaco




**Wire:** more comfortable for the patient under local anaesthetic, particularly when used with a sliding wire speculum or a non-adjustable speculum



**Solid:** holds eyelids securely


## Specula

Non-adjustable – Infant

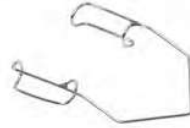

0104090 **SS**

Barraquer speculum for baby with 12mm closed wire blades. 14mm maximum opening, overall length 30mm

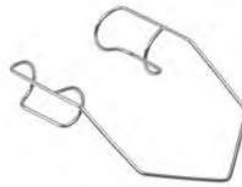

0104075 **SS**

Barraquer speculum for child with 12mm closed wire blades. 15mm maximum opening, overall length 42mm



0104130 **SS**

Barraquer speculum for child with 11mm solid blades. 18mm maximum opening, overall length 40mm


Note: Maximum opening for non-adjustable speculum is measured at rest and may be reduced under lid tension

## Specula

### Non-adjustable – Adult

**0104074** **SS**

Barraquer speculum for adult with 17mm closed wire blades. 16mm maximum opening, overall length 45mm

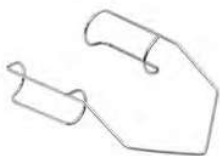

**0104131** **SS**

Barraquer speculum for adult with 15mm solid blades. 18mm maximum opening, overall length 38mm




**0104078** **SS**

Kratz-Barraquer speculum for adult with 13mm open wire blades. 18mm maximum opening, overall length 45mm




**0109401** **Ti**

Kratz-Barraquer speculum for adult with 8mm open wire blades. 18mm maximum opening, overall length 42mm

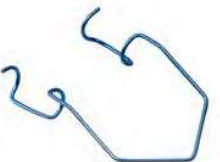

**0104087** **SS**

Kratz-Barraquer temporal incision speculum for adult with 12mm open wire blades, 18mm maximum opening, overall length 45mm

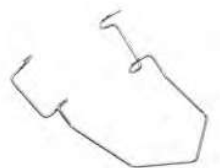

**0104135** **SS**

Sauer speculum with 10mm solid blades. 18mm maximum opening, overall length 16mm



Note: Maximum opening for non-adjustable speculum is measured at rest and may be reduced under lid tension

# Specula

Adjustable – Infant



0104006

Clarke speculum for baby with 8mm open blades.30mm maximum opening, overall length 65mm



0104004

Clarke speculum for child with 12mm open blades.30mm maximum opening, overall length 70mm



0104132

SS

Jaffe lid retractor, small with 10mm open wire blade, overall length 40mm. Comes as a pair.

## Specula

### Adjustable – Adult

**0104085** **SS**

Barraquer sliding speculum for adult with 17mm closed wire blades. 25mm maximum opening, overall length 35mm



**0104083** **SS**

Kratz-Barraquer sliding speculum for adult with 12mm open wire blades. 25mm maximum opening, overall length 32mm



**0104112** **SS**

Pierse sliding speculum for adult with 12mm open wire blades, 25mm maximum opening, overall length 35mm



**0104002**

Clarke speculum for adult with 14mm open blades. 30mm maximum opening, overall length 85mm



**0109410** **Ti**

Lieberman speculum for adult with 15mm open 'K' wire blades. 45mm maximum opening, overall length 84mm



## Specula

Adjustable – Adult



0104142 SS

Lieberman speculum for adult with 15mm open 'K' wire blades. 40mm maximum opening, overall length 82mm\*



0104128 SS

Lieberman speculum for adult with 15mm open 'V' wire blades. 35mm maximum opening, overall length 82mm\*



0104129S SS

Lieberman temporal incision speculum for adult with 15mm open 'K' wire blades. 35mm maximum opening, overall length 78mm\*



0104137 SS

Lieberman speculum for adult with 9mm solid blades. 40mm maximum opening, overall length 82mm



0104136 SS

Lieberman speculum for adult with 15mm solid blades, suitable for both standard and temporal incision 45mm maximum opening, overall length 75mm

\*Aspirating versions are available on request

## Specula

### Adjustable – Adult

**0104141** **SS**

Knapp speculum for adult with 15mm open wire blades. 33mm maximum opening, overall length 84mm



**0104134** **SS**

Lancaster speculum for adult with 15mm solid blades. 38mm maximum opening, overall length 82mm



**0104133** **SS**

Jaffe lid retractor, large with 15mm open wire blade, overall length 40mm. Comes as a pair.

